


Parts of a Microscope

Eye Piece/Ocular Lens

- Magnifies the specimen image
- Where you look through to view the object

Body Tube

- Supports the eyepiece and lenses
- Maintains the proper distance between the eyepiece and objective lenses

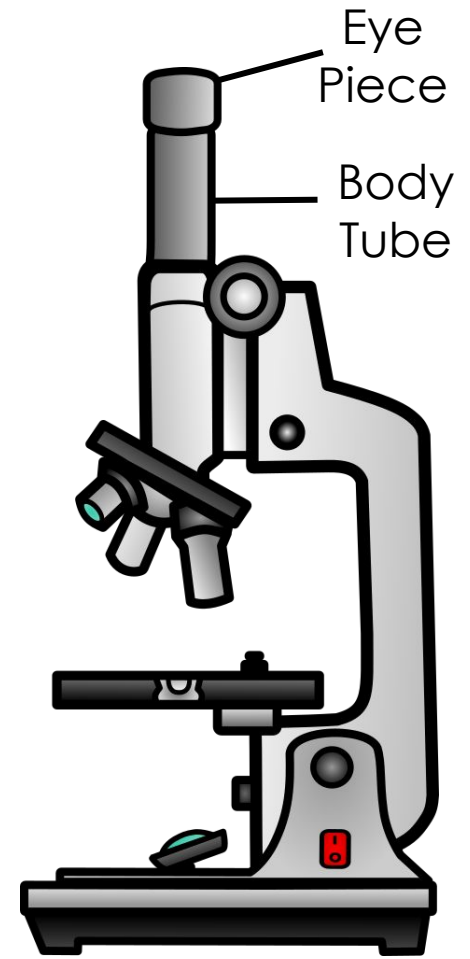

Focusing

1. Always start with the scanning objective
 - Use the coarse knob to focus and then the fine until clear
2. Once it is focused, switch to low power
 - Use the coarse knob to focus and then the fine until clear
3. Switch to high power
 - Only use the fine adjustment to focus

Cleaning and Storing

- Clean only with a soft cloth
- Make sure it is on a flat surface
- Carry with two hands
- Store with the scanning objective in place
- Wrap cords and cover microscope
- Place in the designated area